


RIC hearing loss	profound	7px, 5px Cellion primax	7px, 5px, 3px Pure primax	7bx, 5bx, 3bx Carat binax	Essential Orion 2
	moderate to severe	7px, 5px, 3px Ace primax			
	moderate				
	mild				
BTE hearing loss	profound	7px, 5px, 3px Motion primax (SX, SA and P)	7mi, 3mi Nitro micon	Essential Orion 2	Essential Sirion 2
	moderate to severe				
	moderate				
	mild				
Custom ITE / ITC hearing loss	profound		7mi, 3mi Nitro micon	Essential Orion 2	Essential Sirion 2
	moderate to severe	7px, 5px, 3px Insio primax			
	moderate				
	mild				
Custom CIC / IIC hearing loss	profound		7mi, 3mi Nitro micon	Essential Orion 2	Essential Sirion 2
	moderate to severe	7px, 5px, 3px Insio primax	7px, 5px, 3px Silk		
	moderate				
	mild				

* miniTek only (not compatible with easyTek or the easyTek App) ^ n/a Insio IIC or Orion 2 RIC 10

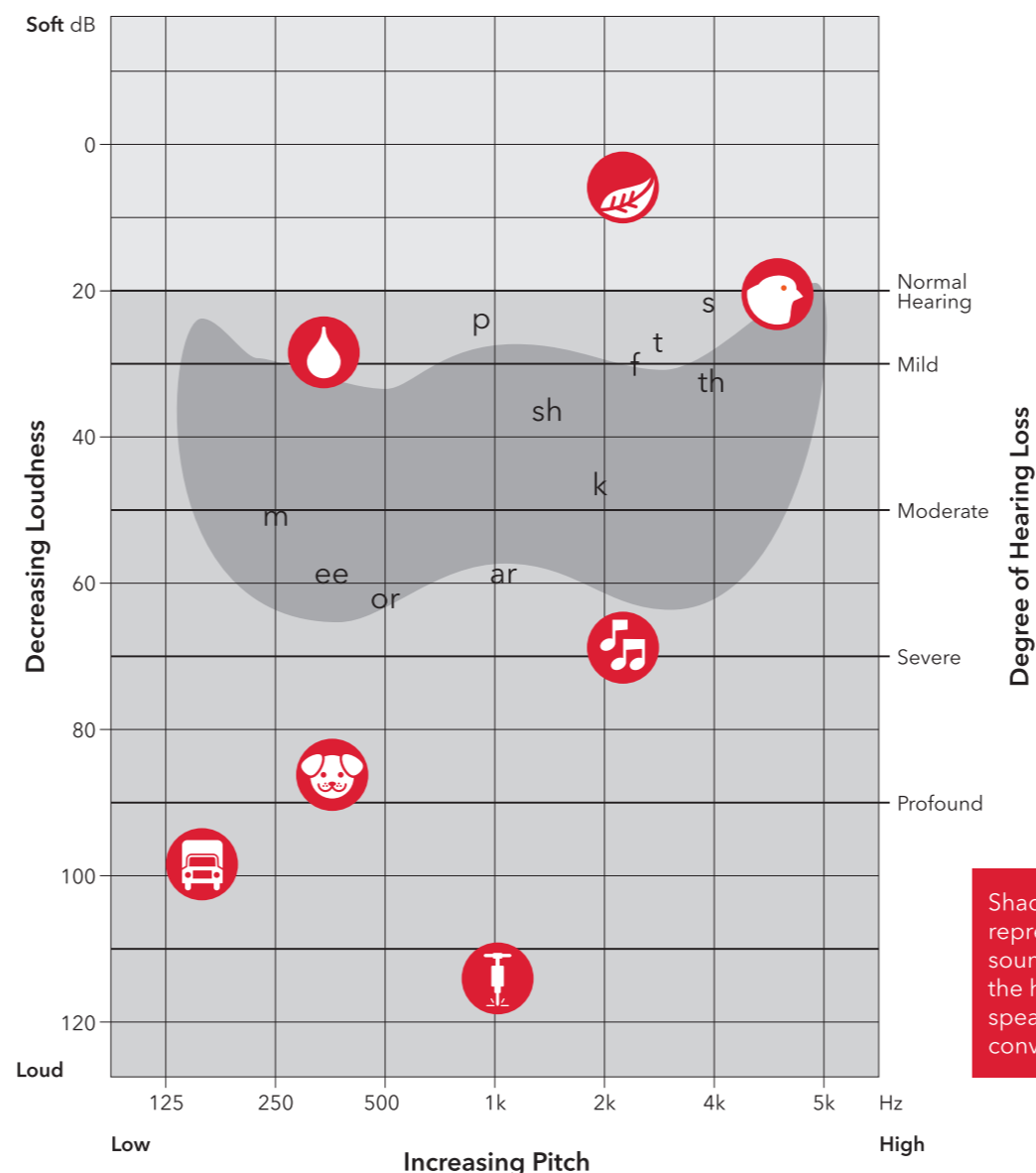

Shaded figure represents all of the sounds that make up the human voice when speaking at normal conversational levels.

RIC	BTE	Custom ITE / ITC	Custom CIC / IIC

*primax CIC features industry first Binaural OneMic directionality

Guide to Signia Hearing Aids



Hearing Systems
SIEMENS

Accessories:



Remote controls: easyPocket and miniPocket.
Our sleek easyPocket™ and miniPocket remote controls are the ideal solution for anyone who wants discreet, convenient hearing aid control without the need for a large remote or smartphone.



Bluetooth compatible: easyTek.
easyTek is a lightweight audio streamer and remote control that quickly connects all e2e wireless™ 3.0 enabled hearing aids to Bluetooth-enabled devices.



Audio Streamer: TV Transmitter
The TV Transmitter wirelessly streams audio from the television to the hearing aids via easyTek. The new TV transmitter is easy to use and supports multiple audio outputs including digital audio.



Companion microphone: VoiceLink.
VoiceLink™ wirelessly streams speech to the hearing aids via easyTek. In difficult listening situations, the VoiceLink™ microphone can be given to speakers to improve speech clarity



Charge Station: Cellion.
The Cellion™ primax inductive charger is used to charge Cellion's™ lithium-ion battery. A single charge provides up to 24 hours use, so your clients never need to worry about their batteries running out on them unexpectedly.



Charge Station: eCharger.
The eCharger recharges and dehumidifies hearing aids overnight, so your clients don't need to worry about changing batteries. Compatible with Pure, Carat and Motion SX/PX models only.



Wireless CROS/BICROS.
The CROS Pure wireless microphone is an energy efficient transmitter for CROS/BICROS Fittings. Delivering exceptional sound quality it is compatible with all primax hearing aid styles using e2e wireless 3.0 technology.



touchControl and easyTek Apps.
The touchControl App offers users discreet control of their hearing aids (e.g. volume/programs) via their smartphone. The easyTek app transforms users' smartphones into an easyTek control center and unlocks advanced hearing aid functionality.

Apps:

Technology to Suit Your Lifestyle

	Essential	Standard	Advanced	Premium	
	Sirion 2	Orion 2	3 primax	5 primax	7 primax
Church / Hall					
Car					
Music					
Outdoor Conversations					
Restaurant / Social Events					
Meetings / Conferences					
Phone Calls					
Watching TV					
1:1 Quiet Conversations					

Key: ■ Advanced Performance ■ Standard Performance □ Basic Performance

Technology Feature Overview

	Essential	Standard	Advanced	Premium	Benefits	
	Sirion 2	Orion 2	5 primax	7 primax		
Church / Hall						
Car						
Music						
Outdoor Conversations						
Restaurant / Social Events						
Meetings / Conferences						
Phone Calls						
Watching TV						
1:1 Quiet Conversations						
	Sirion 2 16 channels/8 handles Auto Directional Mic Speech & Noise Management Feedback Cancellation TouchControl App Frequency Compression Acclimatisation Manager Datalogging Recommended for:	Orion 2 16 channels/8 handles e2e wireless 3.0 easyTek & easyPocket Tinnitus Therapy Signals eWindscreen SoundSmoothing Adaptive Directional Mic Recommended for:	3 primax 24 channels/12 handles Standard SpeechMaster Standard HD Music TwinPhone Standard Narrow Directivity TruEar Learning Tonal Tinnitus Therapy Ocean Waves Therapy Signal Recommended for:	5 primax 32 channels/16 handles Advanced SpeechMaster Advanced Narrow Directivity SpeechFocus Direct. Speech Enhancement Standard Sound Equaliser Spatial Configurator Binaural eWindscreen (man) Recommended for:	7 primax 48 channels/20 handles Premium SpeechMaster Premium HD Music EchoShield Binaural OneMic Directivity Premium Narrow Directivity Spatial SpeechFocus Premium Sound Equaliser Extended Bandwidth Binaural eWindscreen (auto) Recommended for:	Benefits Premium control system that isolates the target speaker from unwanted background noise & speech Three dedicated premium music programs for an outstanding music listening experience Improves listening ease by reducing echo in reverberant situations (e.g. atriums, hallways) Provides frontal directional focus in discreet one microphone (e.g. CIC) custom styles Premium narrow directionality clinically proven to outperform normal hearing* Focuses on speech automatically from any direction Extended fine-tuning of the automatic classification system (up to 6 classes) High-Fidelity sound perception and audibility for frequencies up to 12 kHz Premium automatic binaural solution for wearing comfort in windy situations Advanced control system that isolates the target speaker from unwanted background noise & speech Narrows the directional beam of the directional focus further than standard narrow directionality Focuses on speech automatically from the back or front Reduces unwanted speech noise in cocktail party situations Standard fine-tuning of the automatic classification system (up to 3 classes) Manual control over the width and direction of the directional beam Manual binaural solution for wearing comfort in windy situations Standard control system that isolates the target speaker from unwanted background noise & speech A dedicated music program to enhance music listening Allows you to automatically hear your phone conversations in both ears Narrows the directional beam of the directional focus Emulates the natural localisation function of the outer ear Automatically learns the wearer's amplification preferences Clinically proven patented approach to achieve tinnitus relief without adding masking sound Customisable standard as well as ocean-wave therapy signals for tinnitus relief Synchronises bilateral hearing aids and facilitates communication to the easyTek and easyPocket Wirelessly connects your hearing aids to different Bluetooth®-enabled devices & remote controls Standard customisable therapy signal for tinnitus relief Advanced solution for improved wearing comfort in windy situations Instantly suppresses sudden unpleasant transient noises Multi-channel adaptive suppression of moving sources of background noise Automatically activates frontal directional focus Enhances speech while reducing unwanted background noise Advanced solution for suppressing annoying whistling Free smartphone app for remote control of basic hearing aid functions Allows perception of high-pitched sounds for patients with limited high frequency hearing Adjusts automatically to the wearer's preferred amplification Accurate recording of hearing aid usage e.g. time worn

*Two clinical studies have shown that the narrow directionality feature provides better than normal hearing in certain demanding environments (University of Northern Colorado, 2014; Oldenburg Hörzentrum, 2013): Speech Reception Thresholds (SRT) in cocktail-party situations improved up to 2.9 dB for wearers with mild to moderate hearing loss compared to people with normal hearing.